

Eye Irritation Test: Potency Assessment Method (ET50) for Mild and Ultra-Mild Test Items

Test code: CT-027

Quality Statement

XCELLR8 is accredited by the UK Medicines and Healthcare Products Regulatory Authority (MHRA) for the conduct of *in vitro* safety testing in compliance with Good Laboratory Practice (GLP). This means that we are able to provide our clients with test results that may be used at a regulatory level to demonstrate product safety, where the test is an approved regulatory method.

Overview


Eye irritation is defined as changes in the eye following the application of a test chemical to the surface, which are fully reversible within 21 days of application.¹

The test described here is a non-regulatory method for the potency assessment of mild and ultra-mild chemicals in terms of human eye irritation potential. It provides discrimination between Mild, Milder and Mildest formulations.

The test is based on the depth of injury model or Maurer and Jester, using the rationale that the degree of irritation caused by a test substance correlates with the degree of penetration into cell layers and subsequent impact on cell viability.

The method utilises reconstructed human cornea-like epithelium (RhCE), which in its overall design mimics the biochemical and physiological properties of the corneal epithelium of the human eye. The test item is applied directly to the cornea surface, providing a good model of "real life" exposure. Cell viability is measured by enzymatic conversion of the vital dye MTT into a blue formazan salt that is quantitatively measured after extraction from the skin tissues. Irritation potential is calculated in terms of the "ET50" value: the time taken, in minutes, for the test item to reduce the viability of the skin model to 50%. ET50 values are then used to assign the irritancy classification based on the proven prediction model.

Test System: Reconstructed Human Cornea-like Epithelium

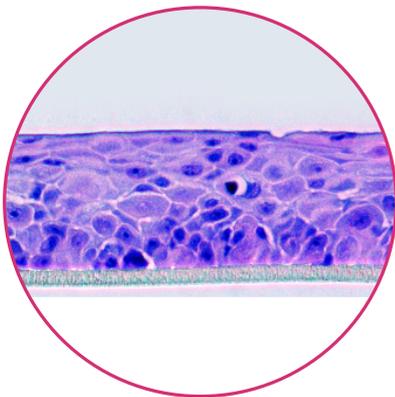

Cross section through reconstructed human cornea-like epithelium



Application of a test item to the model surface




MTT test

Reconstructed human cornea-like epithelium (EpiOcular™) is a corneal model composed of living human cells which have been cultured to form a multi-layered, differentiated corneal epithelium. The levels of differentiation obtained are at the cutting edge of *in vitro* tissue technology. The model consists of highly organized basal cells which progressively flatten out as the apical surface of the tissue is approached, analogous to the normal human *in vivo* corneal epithelium. The profiles of key differentiation markers also mirror those seen *in vivo*. The cells are both metabolically and mitotically active, and release many of the pro-inflammatory agents (cytokines) known to be important in eye irritation and inflammation. EpiOcular™ is grown on special platforms at the air-liquid interface, allowing for direct application of test items in a way that accurately models "real life" eye exposure.

+44 (0)1925 607 134 | info@x-cellr8.com | www.x-cellr8.com

XCELLR8 Ltd. Registered in England and Wales. Registered office: The Innovation Centre, Sci-Tech Daresbury, Keckwick Lane, Daresbury, Cheshire, WA4 4FS, UK. Company number 6489519. VAT number GB 932 3310 59.

 **XCELLR8**
replacing animal testing

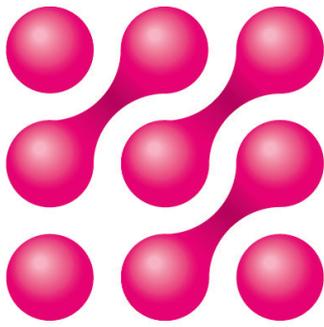

Eye Irritation Test (cont): Potency Assessment Method (ET50) for Mild and Ultra-Mild Test Items

Test code: CT-027

Turnaround Time

4 – 6 weeks

Amount of Sample Required

Minimum 10ml (liquids) / 10g (solids)

Price

Our test prices are dependent on the quantity of test items. Please enquire for a quote using the contact information shown below, or the contact form on our website.

Summary of the Test Method

- EpiOcular™ models are pre-warmed in a cell culture incubator (37°C / 5% CO₂) for 60 minutes or overnight. The culture medium is replaced prior to applying treatment.
- The test item is applied to the surface of the EpiOcular™ models: triplicate models are dosed at the apical surface with 100µl.
- Controls consist of ultrapure water (negative control) and 0.3% Triton X-100 (positive control).
- The dosed skin models are placed into a cell culture incubator for 1 hour, 5 hours and 24 hours (triplicate models per time point).
- Test items and control substances are removed from the skin models surface by washing.
- Following a post-exposure soak in culture medium, the viability of the EpiOcular™ models is assessed by MTT conversion. MTT solution is applied to the surface of the models and placed into a cell culture incubator for 3 hours. The blue formazan metabolite produced by viable cells is then extracted into isopropanol by incubation at room temperature for 2 hours.
- Triplicate samples of the extracted formazan solution are transferred to a microplate and the formazan product is quantified by absorbance spectrophotometry (wavelength 570nm).
- Absorbance readings of the formazan product from EpiOcular™ models incubated with test items are compared with those of negative controls to calculate percentage viability. Absorbance values from all 3 time points are used to calculate ET50 (the time, in minutes, taken to reduce the viability of the EpiOcular model to 50% of the negative control value).
- ET50 values are placed into rank order and / or compared with those of known benchmarks to provide a prediction of the eye irritation potential of the test item.
- A range of acceptance criteria must be satisfied in order for the experimental run to be valid.

Further Downloads

[XCELLR8 Good Laboratory Practice \(GLP\) Compliance Certificate.](#)

References

UN (2009), United Nations Globally Harmonized System of Classification and Labelling of Chemicals (GHS), Third revised edition, UN New York and Geneva.

+44 (0)1925 607 134 | info@x-cellr8.com | www.x-cellr8.com

XCELLR8 Ltd. Registered in England and Wales. Registered office: The Innovation Centre, Sci-Tech Daresbury, Keckwick Lane, Daresbury, Cheshire, WA4 4FS, UK. Company number 6489519. VAT number GB 932 3310 59.

